
INSTRUCTION MANUAL


SMART BODY FAT SCALE A520BT





THANK YOU FOR YOUR PURCHASE.

Please do read the user manual carefully and thoroughly so as to ensure the safe usage of this product, and keep the manual well for further reference in case you have problems.



The Smart Body Fat Scale offers data of your personal weight, body fat percentage, body water percentage, muscle mass; bone mass, BMI (Body Mass Index) and more.

The scale Provides you an accurate measurements of your overall health and fitness.



IMPORTANT SAFETY PRECAUTIONS:

When using digital scale, certain basic safety precautions must be followed, including the following:

1. Scale is a high precision measuring device. Never jump or stomp on the scale or disassemble.
Please handle the scale carefully to avoid breakage.
2. To ensure weighing accuracy, when the scale is first used, please place them on an even surface and wait at least 10 seconds before stepping on the scale. It also applies if the scale has been moved. We suggest you to taking your measurements at the same time every day to ensure the accuracy.
3. Check the device before each use. Do not use the device or an electrode if it is damaged in any way.
The continuous use of a damaged unit may cause injury, improper results, or serious danger.
4. Please check the batteries if the scale malfunction. Change to new batteries if needed.
5. When it is not in use, please remove the batteries and place the device in a dry room and protect it against extreme moisture, heat, lint, dust and direct sunlight. Never place any heavy objects on the equipment.
6. The scale contains plastic materials, if you are Allergic to plastic, do not use it.



For this Bluetooth Scale, there are other safe precautions must be followed, including the following:

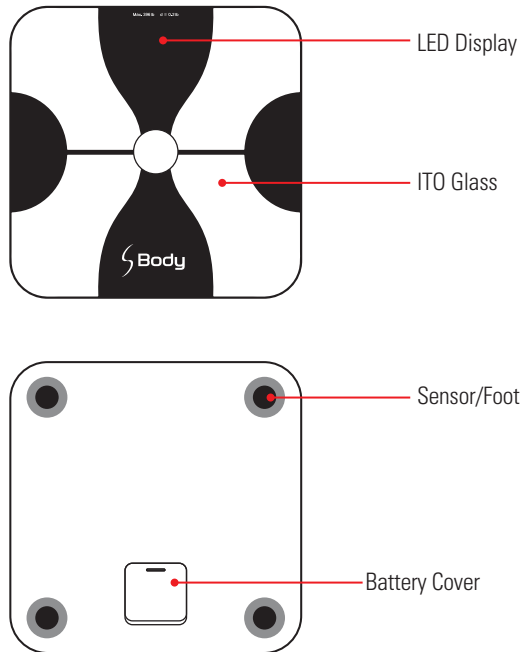
1. Use the Scale while barefoot. The scale cannot measure body fat percentage if you're wearing shoes or socks.
2. Make sure your feet are dry before weighing yourself.
3. **Place the scale on a hard, flat surface during use. Using the scale on a soft surface - like carpet or linoleum - may affect the accuracy of the Scale readings.**
4. The Bluetooth scale should not be used to diagnose or treat any medical condition. You should consult a doctor when undertaking any diet or exercise program. If you are pregnant/pacemaker or other internal device holder, you should not rely the data provided by this product. Please note that all the data are for reference and is not a substitute for advice from a licensed health-care professional.
5. This scale is not waterproof; never immerse the scale in water. To clean the surface, use a damp cloth or a glass cleaner. Do not use soap or other chemicals.



GENERAL INSTRUCTIONS

The Smart Body Fat Scale applies SENS ON technology. It's not necessary to turn on or turn off the scale, just stand on the scale for measurement. It is delicate and trendy, convenience for using.

DEVICE COMPONENTS



INTENDED USE

The Smart Body Fat Scale is intended only for usage at home.

WHAT IS INCLUDED?

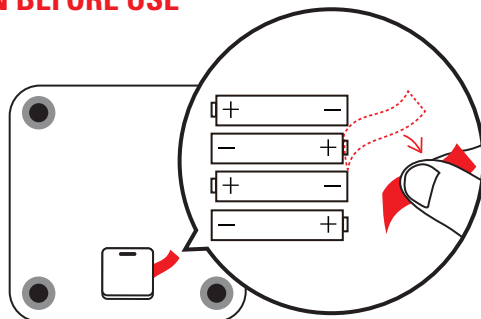
The Smart Body Fat Scale box includes:

- S Body Smart Body Fat Scale
- 4 x AAA Batteries(1.5V each)
- Instruction Manual

LED DISPLAY



PREPARATION BEFORE USE



Note: Remove insulating strip beneath the batteries in the battery holder at the back on the scale.

To SET WEIGHT UNIT

Two units: *lb, kg*

The default unit on scale is decided by the UNIT Selection in APP, please change default unit in APP for daily use. If the scale is with unit systems of Imperial (pounds) and Metric (kilograms), users can select from units lbs to kg.

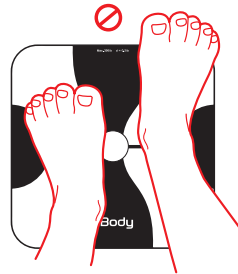
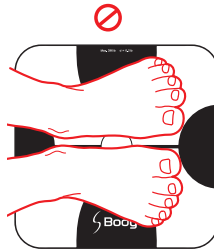
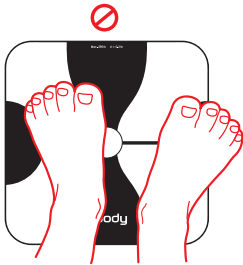
START USING

Set each foot on each side of the center point naturally as show below:

Correct!




Note: To ensure the accuracy please stand on the scale barefoot with the right posture.




OPERATING INSTRUCTIONS


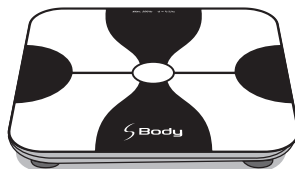
1. Place the scale on hard, flat ground.



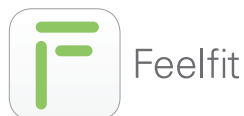
2. Download the App: Search "Feelfit" in the iOS App Store or the Android Application Market and turn on the Bluetooth functions of your device.



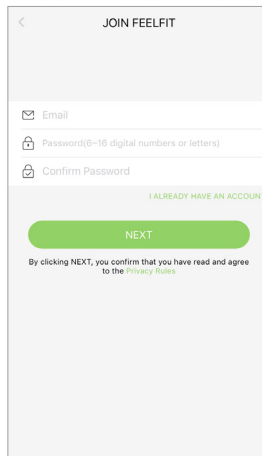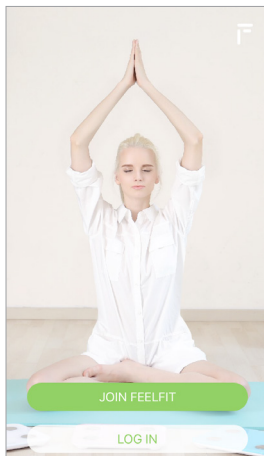3. Register your own ID by email for keeping record. Add personal data during register and confirm OK.



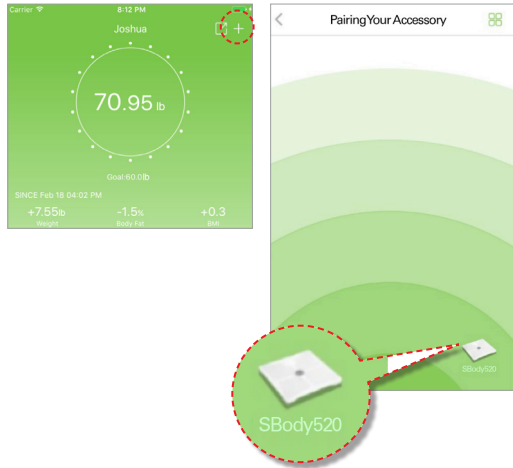4. Press “+” on the measuring page to find device “**SBody520**”.



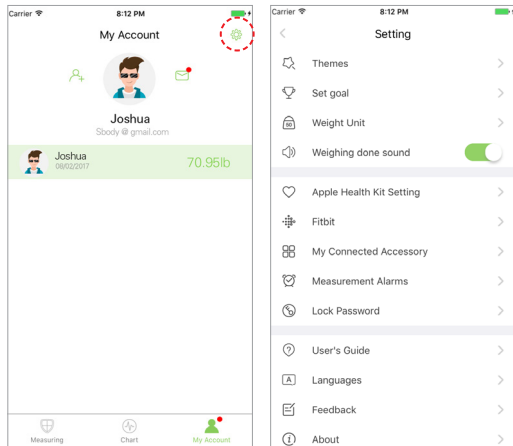
5. Step on the scale with BARE feet and wait for 10-15 seconds, the circle is stable with “ding” sound when you turn on bell ring, it means weighing and analyzing are finished.




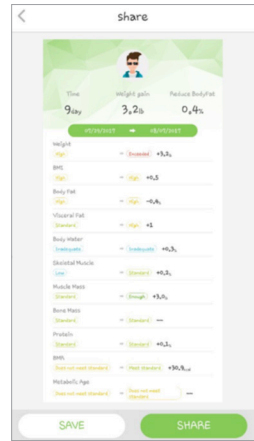
6. Click each data, you can find the analyzing report on body weight, BMI, body fat, water, Skeletal muscle, bonemass, protein, BMR and body age. When you click the “chart” item at bottom, you can see your long-term tracking in chart form.




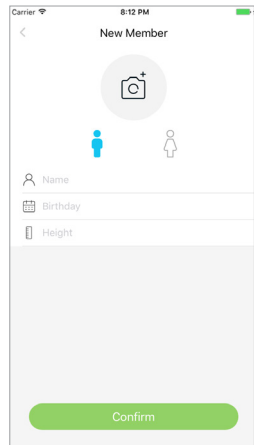
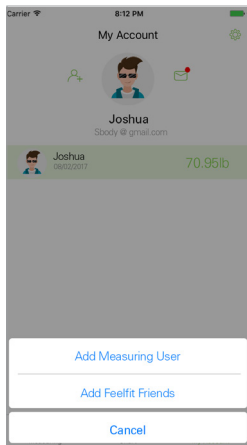
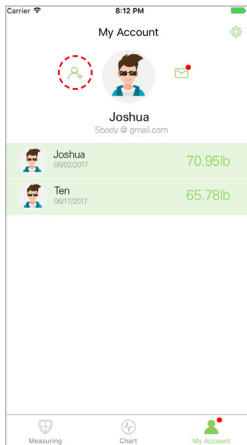
7. When you enter “My Account Setting”, Press “⚙️”, Themes, Set goal, Weight Unit, Weight done sound, languages, Password, and etc. can be set on My Account page. And you can switch the weight unit by pressing setting and entering Weight Unit.



8. Press “”, you can share results of your data to someone.



9. When you enter “My Account”, you can add unlimited user member by press “” and also add new users’ personal data.



SPECIFICATIONS

Product Name	S Body-A520BT
Dimension(W x D x H)	11.96 x 11.96 x 0.86inch / 304 x 304 x 22mm
Net Weight	4.1lb / 1.9kg (Including the battery)
Display	Highlight LED (3.15 x 1.22inch / 80 x 31mm)
Measurement Weight Unit	lb (Pounds) / kg (Kilograms)
Measurement	Body Weight, BMI (Body Mass Index), Body Fat, Water Skeletal Muscle, Bone Mass, Protein, Calorie and Body age
User Number	Unlimited user member
Measurement Range	11lb ~ 396lb / 5kg ~ 180kg
Division	0.2lb / 0.1kg
Accuracy	0.2lb / 0.1kg
Working Environment	32°F ~ 104°F / 0°C ~ 40°C Humidity: ≤ 90% RH
Storage Environment	-4°F ~ 140°F / -20°C ~ 60°C
Power Source	6V (4 x AAA size Batteries)
Auto-ON	SENSE ON Technology
Auto-OFF	After 30 seconds while showing 0.0 About 35 seconds after the weight data is locked

About the Accuracy of This Product

- This product passed strict inspection before delivery and therefore its accuracy is guaranteed the manufacture. Please refer to the above table for the descriptions on accuracy.
- The typical service life of the new and unused batteries for AAA batteries is 4months measurements for the operation time for AAA batteries is 480 times.
- This product is specially designed for Body Fat Analyze as well as measurement. It should NOT be used by anyone during the process of transaction for verification of goods weight.

WARRANTY

warrants its products free of defects in materials and workmanship in normal use for a period of ONE year from the date of retail purchase.

To ensure speedy handing of your issue, please email us for assistance.

Email: contact@visiontechshop.com



VISIONTECH
SHOP
www.VisionTechShop.com